

SHZ-D(III) Two Meters double taps Anti-corrosion pump head

PRODUCT STRUCTURE



Technical Parameter

MODEL	SHZ-D(Ⅲ)Double- gauge/Double-tap	SHZ-D(皿)Four- gauge/Four-tap
Power(W)	180	370
Voltage(V/Hz)	220 50/60	
Flow(L/min)	60	
Lift(m)	8	
Body Material	Anti-Corrosive	
Max. Vacuum(Mpa)	0.098	
Sucking Rate for Single Tap(L/min)	10	
Quantity of Suction Air Hole(Pcs)	2	4
Tank Capacity(L)	15	
Remark	Equipped with 2L, 5L glass reactor or rotary evaporator. The vacuum degree is better when the water temperature is between 4℃ and 10℃	

Tips:110V/60Hz is available


PRODUCT DISPLAY





Point 1 ► All Copper Core Coil Motor

180W High-power Motor High temperature resistance, Long working life

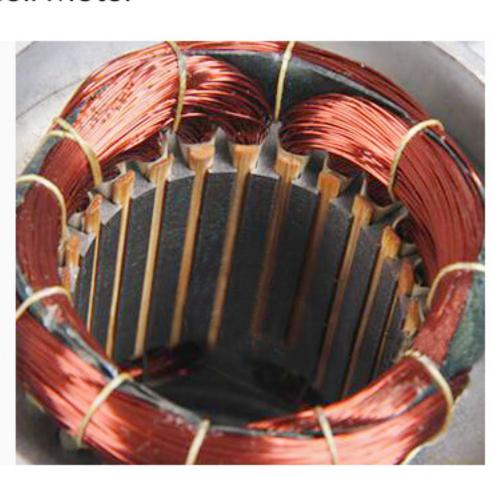

Point 2 > Stainless steel axis of the motor



Anti-corrosion, abrasion resistance and long operating life

Point 3 ► Sink

A molding engineering plastic main part, not easy to deform, durable use with long service time



Point 4 ▶ Independent two meters double taps



Provide vacuum source for two containers

Point 5 Turbine bypass

Six-channel impeller design, great pumping speed and force



DETAILS SHOW

